
ADLER

EUROPE

AD 2176




| | | |
|-------|----------------------------|-------|
| (GB) | user manual | - 5 |
| (DE) | bedienungsanweisung | - 9 |
| (FR) | mode d'emploi | - 14 |
| (ES) | manual de uso | - 19 |
| (PT) | manual de serviço | - 24 |
| (LT) | naudojimo instrukcija | - 29 |
| (LV) | lietošanas instrukcija | - 33 |
| (EST) | kasutusjuhend | - 37 |
| (HU) | felhasználói kézikönyv | - 49 |
| (BS) | upute za rad | - 45 |
| (RO) | Instrucțiunea de deservire | - 41 |
| (CZ) | návod k obsluze | - 63 |
| (RU) | инструкция обслуживания | - 71 |
| (GR) | οδηγίες χρήσεως | - 58 |
| (MK) | упатство за корисникот | - 53 |
| (NL) | handleiding | - 76 |
| (SL) | navodila za uporabo | - 110 |
| (FI) | manuaal ng pagtuturo | - 89 |
| (PL) | instrukcja obsługi | - 120 |
| (IT) | istruzioni operative | - 85 |
| (HR) | upute za uporabu | - 67 |
| (SV) | instruktionsbok | - 97 |
| (DK) | brugsanvisning | - 81 |
| (UA) | інструкція з експлуатації | - 93 |
| (SR) | Корисничко упутство | - 102 |
| (SK) | Používateľská príručka | - 106 |
| (AR) | دليل التعليمات | - 114 |






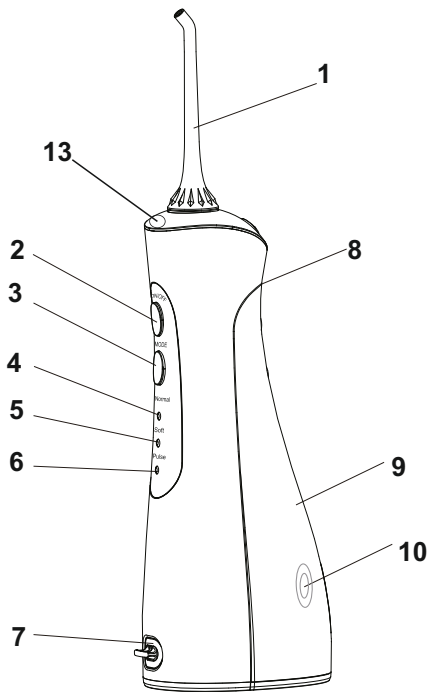
11

12

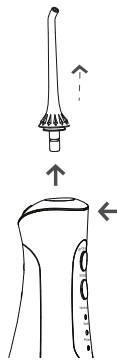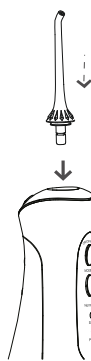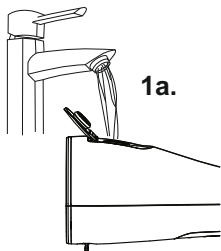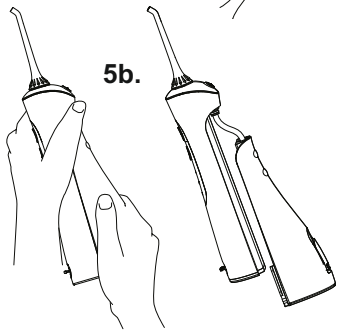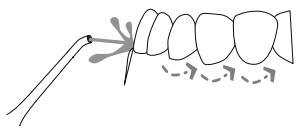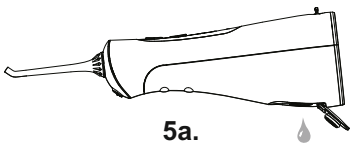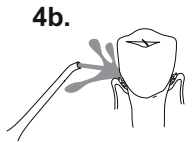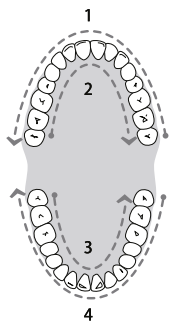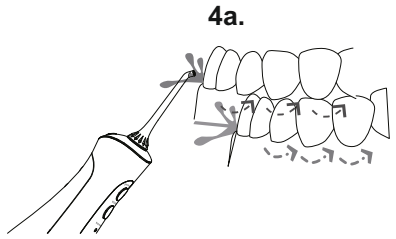SAFETY CONDITIONS IMPORTANT INSTRUCTIONS FOR SAFETY USE PLEASE

PLEASE READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

The warranty conditions vary when the device is used for commercial purposes.

1. Before using the product, please read carefully and always follow the instructions below

instructions. The manufacturer is not liable for any damage resulting from improper use.

2. The product may only be used indoors. Do not use the product for any purpose other than its intended use.

3. The applied voltage is 5V, ~ 1A. For safety reasons, it is not appropriate to connect multiple devices to one power outlet.

4. Be careful when using children. Don't let the kids play product. Do not let children or people unfamiliar with the device use it without

supervision.

5. **WARNING:** This appliance is intended for use by children over 8 years of age and by persons who are assigned to it

diminished physical, sensory or mental capacity, or persons with no experience or

knowledge of the device only under the supervision of the person responsible for the device

safety or if they have been instructed in the safe use of the device and are aware

the dangers of its operation. Children should not play with the device.

Cleaning

and maintenance of the appliance should not be carried out by children unless they are finished

8 years and these activities are carried out under supervision.

6. Always gently pull out the plug after using the product power socket holding the outlet by hand. Never pull on the power cord !!!

7. Never put the power cord, plug or the entire device into water. Never expose the product to the elements such as direct sunlight or rain, etc.

8. Periodically check the condition of the charging cable. If the charging cable is damaged, the product should be referred to a professional

- service center who will be replaced to avoid any dangerous situations.
9. Never use the product with a damaged power cord or if it has been dropped or damaged in any way or if it does not function properly. Do not try to repair a damaged product alone, as this could lead to electric shock. Always have a damaged device replaced by a service center for repair. All repairs can only be done by authorized service technicians. Incorrectly performed repairs may cause dangerous situations for the user.
 10. Never put the product on hot or warm surfaces or kitchen appliances, such as an electric oven or gas burner.
 11. Never use the product near flammable materials.
 12. Do not let the cord hang over the edge of the meter or touch hot surfaces.
 13. Never leave the product connected to the power supply unattended. Even when use is interrupted for a short time, disconnected from the mains, disconnect the power.
 14. For additional protection, it is advisable to install a residual current device (RCD) in the electrical circuit with a rated residual current not exceeding 30 mA. In this regard, please consult a specialist electrician. Contact a professional electrician for this
 15. Install the device in such a way that it cannot fall into water or fall to the floor. The plug should always be available.
 16. Periodically check that the cord is not damaged. If the cord is damaged, take the charging unit to a service center.
 17. A damaged or non-functional unit should no longer be used.
 18. Do not modify or repair the product. It may cause fire, electric shock, or injury.
 19. Cleaning and user maintenance shall not be made by children.
 20. ATTENTION children from 3 to 8 years of age may use the device safely only under the supervision of adults, adolescents from 8 to 14 years of age can safely use the device under the supervision of adults or after having been instructed on the safe use of the device and understand the related dangers, people with reduced physical, sensory or mental abilities or lack of experience and knowledge can use the product after having been given proper supervision or instruction on

how to use the device, and understand the related dangers.

21. Children must not play with the device.

22. The device should only be used with liquids defined in the manufacturer's instructions.

23. Use this product only for its intended use described in this manual. Do not use attachments or chargers that are not recommended by the manufacturer.

24. If the product is dropped, the brush / irrigator head must be replaced before the next use, even if no damage is visible.

25. Small parts can be a choking hazard.

26. Do not place the device in water or liquid, or store it in a place where it may fall or be pulled into a tub or sink.

27. Do not open or disassemble the device. For battery recycling, please dispose of the entire device according to local environmental regulations.

29. Opening the device will destroy the device and will void the warranty.

30. Always hold the power plug in place of the cord when unplugging.

31. Do not touch the power plug with wet hands. It may cause electric shock.

32. If you are undergoing any oral hygiene treatment, consult your dentist before use.

33. This device is a personal care device and is not intended for multi-patient use in a dentist's office or institution.

34. Do not put water under your tongue, ear, nose, or other delicate areas. This product can build up pressures that can cause serious damage in these areas.

35. Do not aim water at the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, can be present in some tap water or non-chlorinated well water, and can be fatal if directed into these areas.

Device description:

- | | | |
|----------------------------------|------------------------|---------------------|
| 1. Cleaning tip | 2. Power switch | 3. Mode selection |
| 4. Normal mode | 5. Delicate mode | 6. Pulse mode |
| 7. charging socket and split pin | 8. Air supply opening | 9. Water container |
| 10. Flap with water inlet | 11. Spare cleaning tip | 12. Charging cable. |
| 13. Cleaning tip release button | | |

Preparing the device for work:

Before turning on the device, charge it by removing the pin from the charging socket (7). Insert the

appropriate end of the charging cable (12) into the charging socket (7). Do not use the device while it is charging. After completing the charging, insert the pin into the charging socket, making sure that there is no water or other liquid in the socket that may cause a short circuit.

Then insert the cleaning tip into the hole at the top of the device with the wider part facing the hole (2a). When you hear a click, the device is properly installed.

After installing the tip, open the water inlet flap (10), be careful to do it gently and not to damage the flap. Fill the device with water and then close the flap so that the water does not escape.

Detaching the Cleaning Tip:

To change the cleaning tip, press the release button (13) and pull the tip straight up away from the device (2b).

WARNING!

The cleaning tip can be rotated 360 degrees. Grabbing the base of the tip, you can slowly turn it left and right.

Usage:

The water jet should be directed only at the teeth. At a slight angle. Do not put too much water directly on the gums and tongue as this can damage the mouth and cause bleeding. Please adjust the operating mode to your personal preferences by selecting the appropriate mode by pressing the mode selection button (3). The preferred method of use is shown in the illustrations included in this manual (4a, 4b).

Maintenance:

After finishing work, open the water inlet flap and empty the tank (5a). Then detach the tank gently. To do this, grab the front of the device, taking care not to turn on the device, and then pulling down the bottom of the water container, slowly removing the pump hose from the water container.

!! ATTENTION !!

Do not pull too hard on the water tank. Movement should be slow so as not to damage the water hose.

Technical data:

Power Source: 5V \approx 1A

Capacity : 150 ml

Pressure : 30-100 psi / 2.0-6.9 bar

 To protect your environment: please separate carton boxes and plastic bags and dispose them in corresponding waste bins. Used appliance should be delivered to the dedicated collecting points due to hazardous components, which may affect the environment. Do not dispose this appliance in the common waste bin.

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